# Patch iontophorèse articulation

## *Patch iontophoresis*

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*1 box with 1 sachet containing 2 patches of 50 x 70mm.*

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*The patch ARTISOULAGE uses the iontophoresis technique for the relief of muscle and joint pain. In use since the 19th century, iontophoresis, also called ion transfer, is a method that involves the use of a direct current to deliver active through the skin to obtain a greater effect in relief chronic muscle and joint pain.*

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*In order to enhance skin functionality recovery reducing pain several active ingredients have been added, particularly selected to provide the product with further properties.*

*HARPAGOPHYTUM: Has anti-inflammatory and analgesic properties.*

*ARNICA: With anti-inflammatory and analgesic properties, helps regenerate tissue damaged due to contusions or lesions by improving circulation in the area.*

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*Apply a patch on the painful area, allow to act for 8-12 h, remove and replace by a new one.*

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*The measurement of micro-currents and magnetic field generated in vitro are evidenced by the study conducted by Abich Clinical Study, certified company.*

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*The patch for chronic muscle and joint pain is a patch whose main action is providing protection against the external environment and relieving localised pain in joints and the musculoskeletal system.*

*The patch creates a protective barrier against the external environment due to the specific characteristics of its fabric. Its elastic adhesive characteristics provide restriction on movement in the affected area to help prevent further trauma. The protection patch encourages healing.*

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*The ingredients contained within the adhesive matrix support the protective device's main action by providing moisturising, emollient and antioxidant action to the area being treated, to make it easier to remove the patch and avoid unpleasant irritation.*

*- Harpagophytum: Anti-inflammatory activity. Used in natural medicine for the symptomatic treatment of chronic osteoarthritis, back pain and tendinitis. Different clinical studies show the effectiveness of the active ingredient in the treatment of pain associated with arthritis and other musculoskeletal conditions of rheumatic or non-rheumatic origin.*

*- Lavender: Anti-rheumatic activity and anti-inflammatory power. Its soothing and softening action on the skin makes the patch's contact with the skin more comfortable.*

*- Bitter orange: In addition to its well known pleasant scent (orange blossom and neroli), it has a slight sedative capacity, which decreases pain perception.*

### Statut

Dispositif médical classe I.

*Medical device class I.*

# Patch Hydrogel articulations

## Hydrogel Patch

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1 sachet contenant 4 patches de 80 x 90 mm.

*1 printed sachet containing 4 patches of 80 x 90 mm.*

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*Thanks to the special characteristics of the patch's non-woven fabric film, on placing the patch on the skin it exerts a slight immobilising action on it and, therefore, on the associated tissue, thereby impeding movements causing pain. This effect is provided by the special characteristics of the fabric from which the patch is made, which has a high resistance to both horizontal as well as vertical stretching.*

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*Menthol provides a deep cooling sensation to the skin.*

*Eucalyptus is a plant used as a pectoral decongestant in traditional medicine due to its soothing properties. Provides cooling and a degree of muscle relaxation.*

*Arnica is a plant used to repair bruises. Helps regenerate tissue damaged due to contusions by improving circulation in the area.*

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*Ensure that the area where the patch will be applied is clean and dry, to ensure perfect adhesion. The patch should remain in place for 1 hour and reapply in the same area if symptoms persist.*

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*1 foil sachet containing 6 patches (60mm diameter).*

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*Localized excessive sweating or hyperhidrosis has long been neglected.*

*The electrodynamic patch consists of a tissue with integrated photovoltaic cells that use iontophoresis as a mechanism to treat hyperhidrosis. The micro-currents generated by the patch in contact with the skin moderates the activity of the sweat cells in the affected area, helping to reduce excessive sweating.*

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***The iontophoretic patch is an improvement on conventional patches due to its more affective dual action :***

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*Apply 1 patch on the affected area (feet, hands or underarms). The area should be clean, dry and hairless. Leave one night (8 hours) and remove in the morning. Repeat the application the next night.*

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
*It has been shown for many years that the application of weak currents would be the most effective treatment for hyperhidrosis. Many articles have been written which testify to it. (More information on request)*

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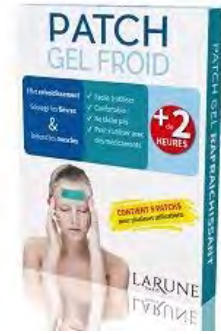*Box with 4 sachets with 1 hydrogel patch each one. TOTAL: 4 hydrogel patches 40 x 110mm*

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*This patch is an advanced hydrogel strip, specially formulated to provide instant cooling relief for forehead and neck, to help to relief ease tight muscles, headaches and migraines, as well as fever.*

*Aid Fresh contains natural menthol which helps to relax tense muscles around neck and head. The gel gives instant relief and provides a cool calming effect the skin. Also is indicated in case of sunstrokes to calm and refresh the affected area.*

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*The refreshing effect of natural menthol achieves:*

- *Quick relief of fevers.*
- *Relief of headaches and migraines.*
- *Helps to relax tense muscles.*
- *Refresh and calms in cases of sunstrokes.*

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*Apply the patch on the forehead and leave on for 1-2 hours.*

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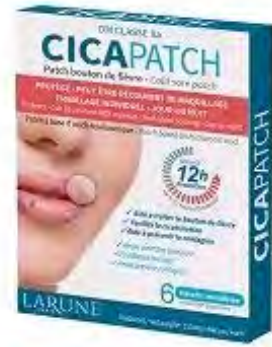
*1 box containing 6 sachets containing 1 patch each one with a leaflet.  
Patch size: 15mm diameter*

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*The cold sore patch with hyaluronic acid is a medical device promoting a faster and better healing of herpes. This transparent patch helps treating Herpes simplex virus facilitating the physiological recovery thanks to their occlusive and mechanical barrier action, mainly intended to keep under control the micro-environment of the infected part.*

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*The cold sore patch contains a salt of hyaluronic acid, a natural component of the skin, necessary for repairing the dermis.*

*The hyaluronic acid is well-known by their moisturizing and soothing effect, which promotes to increase the natural healing process. Add this substance helps the patch to accelerate the vanishing of the cold sore.*

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*Wash hands before the application and be sure that the cold sore area is clean, dry and free of cream in order to grant a good adhesion. Take the patch out from its pouch and put it on. On average each patch can remain in place on dry intact skin for up to 8-12 hours.*

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*1 printed foil sachet containing 1 sheet with 10 films 15 mm in diameter.*

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*The mouth ulcers patch is an innovative health care product (Medical Device) that represents an evolution in the use of low molecular weight hyaluronic acid (hyaluronan), a polysaccharide widely present in body tissues, in which it performs many different fundamental biological functions and through which, thanks to its film-forming and protective properties, it achieves a safer, natural cell regeneration process in case of changes in the oral mucosa. Its mechanism of action is to recover the integrity of the epithelial tissue of the oral cavity (oral mucosa), by using very low viscosity, low molecular weight hyaluronic acid which performs protective activity thanks to its film-forming capacity.*

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*The role of hyaluronic acid in the tissue repair process is well documented in scientific literature. The innovation of this special type of hyaluronic acid is that, thanks to its physical characteristics, it has demonstrated a higher regenerative capacity than the other types of hyaluronic acid and, therefore, provides greater efficacy in process to repair damaged or wounded mucosa.*

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*The patch must be applied at least 3 times a day, in order to maintain the protective film.*

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*1 sachet containing 12 patches, 19 mm in diameter each.  
Sachet size: 100 mm x 75 mm*

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*The warts patch is indicated for the treatment of common warts usually appearing on the hands and feet.  
It is also effective in treating plantar warts.  
Not aggressive for the surrounding skin, it is also suitable for use in children.*

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*The patch contains no active ingredients but uses two mechanisms of action:  
- a mechanical barrier and a perfect occlusion of the skin affected by the wart.  
- a waterproof and transparent patch, which allows at the same time a good maintenance as well as the control of the progression.*

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*Wash your hands before applying. Apply the patch on a clean and dry skin. The occlusive action needs to last for at least 6 days during which the patch should be replaced every 12 to 24 hours.*

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*1 sachet containing 24 patches of 1.5 cm diameter each.  
Size of the bag: 130 mm x 90 mm*

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*This patch helps to treat insect bites. The protection of the patch isolates the skin affected by the insect bite of the external agents, so that by avoiding scratching the bite, the skin does not over-infect. The patch contains active ingredients moisturizing, soothing and refreshing.*

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*The following active ingredients:*

- *Liquorice: soothing and relaxing, it promotes and speeds up the healing process of the skin.*
- *Bisabolol: an active ingredient found in the essential oil of chamomile . It has anti-inflammatory properties to protect the skin.*
- *Menthol: an active ingredient extracted from mint. It has well-known cooling and soothing properties (it also has some local analgesic actions). Antiseptic, it reduces bacterial growth.*

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*Apply directly to the bite and leave between 4 and 8 hours.*

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*1 box with a sachet containing 2 patches, 10 x 20cm.  
Sachet size: 14cm x 24cm*

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*Iontophoretic Patch Scar Reducer is an electrodynamic patch with alternating voltaic cells that uses the principle of iontophoresis to enhance significantly the absorption of active molecules allowing to achieve a greater effect in scar reduction.*

*Iontophoretic Patch Scar Reducer represents the evolution of conventional patches and is much more effective thanks to its double action:*

1. *The electromagnetic field created by the patch generates an ion channel which allows a more efficient absorption of the active ingredients migrating through the skin.*
2. *The electro-ionic action creates a "cell massage" enhancing oxygen exchange between cells, revitalising connective tissue. Therefore, Iontophoretic Patch Scar Reducer provides stimuli that increase collagen production, thus enhancing skin elasticity and health, which leads to speeding up the healing process and reducing these unwanted skin marks.*

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*Iontophoretic Patch Scar Reducer is based on iontophoresis to take the active ingredients more effectively to the site of action :*

*- Hyaluronic acid is a natural skin component, necessary for repairing the dermis, and its deficiency is related to an abnormal arrangement of fibroblasts and loss of hydration, key factors in scar formation. Hyaluronic acid helps us keep the optimum moisturising conditions of the skin to be repaired.*

*· Extractum cepae (onion extract): extract that helps improve collagen production at the sites where it is applied, improving markedly the appearance of scars.*